

# PerfectCut®

## Advanced Tool Set For Perfect Aortic Incisions

### PerfectCut Aortic Lancet

Consistent, reliable, repeatable crisscross incisions

For aortotomies such as at the coronary artery proximal anastomotic and arterial cannulation sites

Lancet blade features retractable safety shield

### PerfectCut Rotating Bullet-Tip Punch

Punch utilizes the lancet's uniform incision to reduce the likelihood of lateral side notches caused by inconsistent scalpel blade incisions

Advanced "Bullet-Tip" designed for straight insertions – no "button-hole" technique necessary

Rotating cutting action requires less force and eliminates shearing of delicate tissue

Self-aligning cutter mechanism minimizes mechanical jamming

REF	SIZE	REF	SIZE
RCK40	+ ● 4.0 mm	RCK48	+ ● 4.8 mm
RCK45	+ ● 4.5 mm	RCK50	+ ● 5.0 mm

Kits sold 6 units/case



LANCET 8 in (20 cm),  
10 x 3.3 mm blade  
Retractable blade shield



PUNCH 8 in (20 cm),  
Bullet-tip, long  
Rotating cutting action

### Step One:

Aortic lancet  
creates precise  
crisscross incisions


①

### Step Two:

Bullet-Tip Punch  
inserts easily into  
incisions

②

Solid, durable  
construction  
results in reliable  
performance

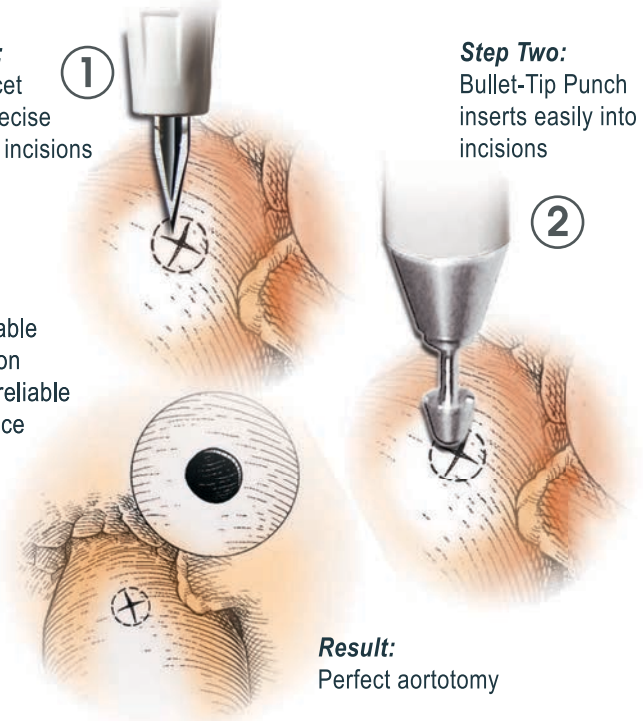

**Result:**  
Perfect aortotomy

 **QUEST Medical, Inc.** [www.questmedical.com](http://www.questmedical.com)

© 2014 Quest Medical, Inc. PerfectCut is a trademark of Quest Medical, Inc. 2010-8

[custserv@questmedical.com](mailto:custserv@questmedical.com)

**R<sub>x</sub> ONLY**

+800.627.0226 (option 1) | +972.390.9800 (v) | +972.390.7173 (f)