

MPS[®] 2

16-CORE HEAT EXCHANGER DELIVERY SETS (5161102)

New larger core for improved cooling. The MPS[®] 16-Core Heat Exchanger provides a larger surface area for faster and more effective cooling. The new 16-Core Delivery Set with Arrest Agent and Additive Cassettes provides another option for your surgical staff looking for colder delivery temperatures¹.

16-Core Delivery Set Includes:


- Crystalloid line
- 16 convolution heat exchanger with new 3rd port for easy access
- Vent line
- 50 ml arrest and additive cassettes
- Blood inlet inner diameter: 6.32 mm
- Blood inlet length: 72 in (183 cm)
- Extension line length 120 in (305 cm)

Capacities of Heat Exchangers

	10-Core	16-Core	Difference
Surface Area	66 in ²	110 in ²	67%
Priming Volume ¹	46 ml	66 ml	20 ml

To retrofit an existing MPS 2 console for the new 16-core heat exchanger, please order an MPS 2 Locking Block Conversion Kit (Cat No. 5001032). It's easy to install and will allow the MPS 2 Console to accept both the 10-Core and 16-Core Heat Exchanger.

For ordering information or questions, please contact your sales rep or customer service. Customer Service: **1.800.627.0226 (option 1)** or custserv@questmedical.com

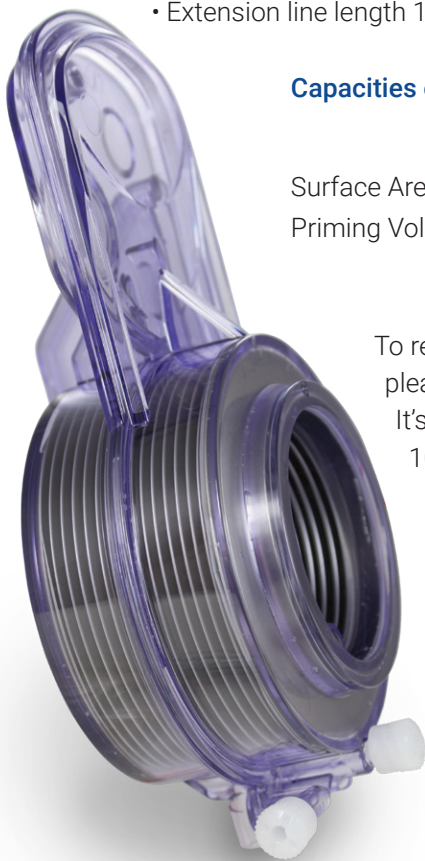

 **QUEST Medical, Inc.**



Rx ONLY

custserv@questmedical.com | 800.627.0226 (option 1)
+1 972.390.9800 | fax: +1 972.390.7173 | www.questmedical.com

¹ Data on file

© 2018 Quest Medical, Inc. MPS is a registered trademark of Quest Medical, Inc.
2018-3