

PerfectCut®

Advanced Tool Set For Perfect Aortic Incisions

PerfectCut Aortic Lancet

Consistent, reliable, repeatable crisscross incisions

For aortotomies such as at the coronary artery proximal anastomotic and arterial cannulation sites

Lancet blade features retractable safety shield

PerfectCut Rotating Bullet-Tip Punch

Punch utilizes the lancet's uniform incision to reduce the likelihood of lateral side notches caused by inconsistent scalpel blade incisions

Advanced "Bullet-Tip" designed for straight insertions – no "button-hole" technique necessary

Rotating cutting action requires less force and eliminates shearing of delicate tissue

Self-aligning cutter mechanism minimizes mechanical jamming

REF	SIZE	REF	SIZE
RCK40	+ ● 4.0 mm	RCK48	+ ● 4.8 mm
RCK45	+ ● 4.5 mm	RCK50	+ ● 5.0 mm

Kits sold 6 units/case



LANCET 8 in (20 cm),
10 x 3.3 mm blade
Retractable blade shield



PUNCH 8 in (20 cm),
Bullet-tip, long
Rotating cutting action

Step One:

Aortic lancet
creates precise
crisscross incisions


①

Step Two:

Bullet-Tip Punch
inserts easily into
incisions

②

Solid, durable
construction
results in reliable
performance

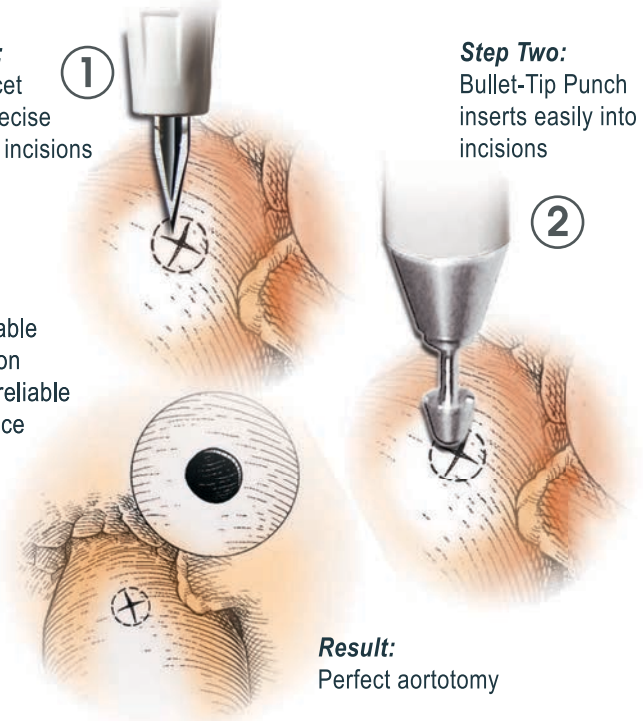

Result:
Perfect aortotomy

 **QUEST Medical, Inc.** www.questmedical.com

© 2014 Quest Medical, Inc. PerfectCut is a trademark of Quest Medical, Inc. 2010-8

custserv@questmedical.com

+800.627.0226 (option 1) | +972.390.9800 (v) | +972.390.7173 (f)

R_x ONLY